

How to use the Bemfola pre-filled pen




You have been prescribed Bemfola and this instruction leaflet will guide you how to use the pen.

Bemfola 75 IU/0.125 ml Bemfola 150 IU/0.25 ml Bemfola 225 IU/0.375 ml Bemfola 300 IU/0.50 ml Bemfola 450 IU/0.75 ml Solution for injection in a pre-filled pen.



Things to
remember:Wash your hands
Give yourself the injection around the same time every day
Take your pen out of the fridge 5–10 minutes before using it

The different parts of your pen

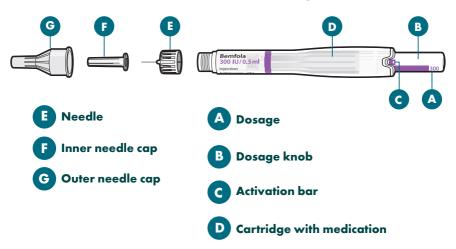

Getting your pen ready

Remove the peel tab from the injection needle. Hold the pen by its sides and attach the needle by screwing it onto the pen.

Caution: Do not push the dosage knob in while attaching the needle.

Correct position of the needle

Incorrect position of the needle






Remove the outer needle cap. Keep it for later. You will need it after the injection. Remove the inner needle cap.






Hold the pen so the needle is pointing upright.¹ Push the dosage knob in until it stops. The activation bar disappears and some liquid should splash out – this is normal.² If no liquid splashes out, do not use the pen, take a new one.

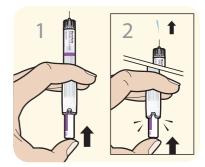

Setting the dose

Turn the dosage knob until your prescribed dose is in line with the display window. The pen is now ready for injection.

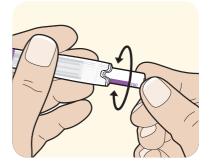

Giving yourself the injection


Immediatly inject the liquid (normally in your tummy) according to the instructions from your doctor or nurse. Lightly pinch the skin of the injection area and insert the needle. Caution! Do not push the dosage knob while inserting the needle!

Push in the dosage knob slowly and continuously until it stops and the dose bar has disappeared. Do not remove the needle immediately. Wait for 5 seconds before pulling it out.


After the injection: Replace the outer needle cap and put the used needle in a closed container (bottle or milkbox) and throw it away in the household waste.

The rest of the used equipment is also thrown in the household waste.



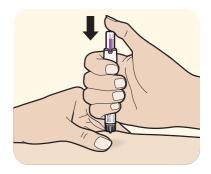

Keep refrigerated ($2^{\circ}C - 8^{\circ}C$). Do not freeze. Can be taken out of the refrigerator and kept at room temperature (max $25^{\circ}C$) for 90 days. Do not use this medicine after the expiry date.

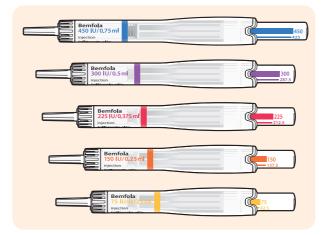

Bemfola 75 IU 0,125 ml	Bemfola 150 IU 0,25 ml	Bemfola 225 IU 0,375 ml	Bemfola 300 IU <i>0,5 ml</i>	Bemfola 450 IU <i>0,75 ml</i>
75 IU	150 IU	225 IU	300 IU	450 IU
62,5 IU	137,5 IU	212,5 IU	287,5 IU	425 IU
50 IU	125 IU	200 IU	275 IU	412,5 IU
37,5 IU	112,5 IU	187,5 IU	262,5 IU	400 IU
	100 IU	175 IU	250 IU	375 IU
	87,5 IU	162,5 IU	237,5 IU	350 IU
	75 IU	150 IU	225 IU	337,5 IU
				325 IU
				300 IU

Bemfola (follitropin alfa), solution for injection in pre-filled pen. Prescription medicine. What Bemfola is used for: Adult women: 1) to help release an egg from the ovary (ovulation) in women that cannot ovulate and that did not respond to treatment with a medicine called "clomiphene citrate", 2) together with another medicine called "lutropin alfa" ("luteinising hormone" or LH) to help release an egg from the ovary (ovulation) in women undergoing very little gonadatropins (FSH and LH), 3) to help develop several follicles (each containing an egg) in women undergoing assisted reproductive technology procedures (procedures that may help you to become pregnant) such as "in vitro fertilisation", "gamete intra-fallopian transfer" or "zygote intra-fallopian transfer". Adult men: together with another medicine called "human chorionic gonadotropin" to help produce sperm in men that are infertile due to a low level of certain hormones. Do not use Bemfola: if you are allergic to follicle stimulating hormone or any of the other ingredients of this medicine, or have a tumour in your hypothalamus or pituitary gland. Do not use Bemfola if you are a woman with large ovaries or sacs of fluids within the ovaries of unknown origin; unexplained vaginal bleeding; cancer in your ovaries, womb or breasts; condition that usually makes normal pregnancy impossible, such as ovarian failure (early menopause), or malformed reproductive organs. Men with damaged testicles that cannot be healed should not use Bemfola. Please read the package leaftet before using Bemfola. This text is based on the package leaftet dated 2024-01-09. Marketing Authorisation Holder and Manufacturer: Gedeon Richter Plc, Gyömröi út 19-21, HU-1103 Budapest, Hungary. Information provided by Gedeon Richter Nordics AB, Barnhusgatan 22, 111 23 Stockholm. *Email: medinfo.se@gedeonrichter.eu*. Telephone: +46 (0)8 611 24 22.

🕞 Gedeon Richter